


The Macy Catheter Bedside Care Kit

Quick guide for insertion, medication administration, and troubleshooting.

Macy Catheter Insertion

1. Using gloves, lubricate catheter tip.
2. Insert to placement marker line between arrows.
3. Using a luer syringe, inflate the balloon with 15 ml of water via the balloon inflation port.
4. Place in between legs and secure to thigh or abdomen.

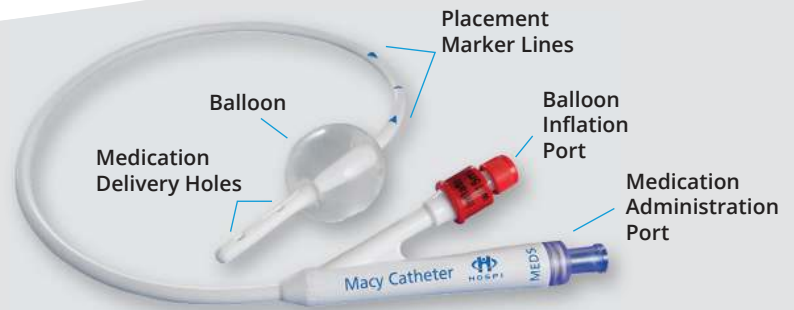

The Macy Catheter®

Medication Administration* Using the Liquipill System

1. Grind Tablets

- Place tablet in slotted opening.
- Turn grinder back and forth until it moves smoothly.

2. Add Water

- Remove grinder cap from medication receptacle and replace with syringe adapter cap.
- Draw up 10 ml of tap water with enteral syringe.
- Insert 10 ml of water into the medication receptacle via syringe adapter cap.

3. Shake/Swirl

- Shake/swirl for 10 seconds.

4. Draw/Inject

- Draw medication solution into 10 ml enteral syringe.
- Inject into medication administration port over 3-5 seconds.

5. Flush

- Flush port by injecting 3 ml of water with 3 ml enteral syringe over less than a second.



Liquipill System

Care and Troubleshooting

1. The Macy Catheter is approved to remain in the rectum up to 28 days.
2. Per your facility's protocol, routine perineal care may be performed with the Macy Catheter in place.
3. If leakage occurs from the Medication Administration Port following medication administration, using the 3 ml enteral syringe, quickly inject 3 ml of water over less than a second.

For more information, visit MacyCatheter.com.
For 24/7 Clinical Support, call 888.509.3732.

*Determination of medications, dosage, and route of delivery is the responsibility of the prescribing clinician.

Patient-centric products that improve quality of life. **HOSPI**



©2018 Hospi Corporation. All rights reserved.